

Hospital	_____
Department	_____
Clinician	_____
Contact information	_____

Provox® LaryTube™

When and How to Use

The Provox LaryTube is a soft silicone cannula designed to hold the Provox Heat and Moisture Exchanger (HME) and Provox FreeHands HME®. It is placed in the stoma to prevent the stoma from shrinking. It can also be used instead of an adhesive when the skin is irritated or during radiotherapy when it is not possible to use an adhesive.

The Provox LaryTube is available in three versions:


- **Fenestrated** – For voice prosthesis users. Designed with holes so that air can go through the voice prosthesis. Can be attached with a Provox TubeHolder or Provox LaryClip™.
- **Standard** – Made for use with or without a voice prosthesis. For voice prosthesis users the LaryTube can be custom fenestrated by your healthcare professional to align with your prosthesis for maximum air flow. Custom punctures can be made using the Provox Fenestration Punch.
- **Blue Ring** – Can be only attached with a Provox Adhesive. Custom punctures can be made using the Provox Fenestration Punch.

1. Inserting the Provox LaryTube:

- Lightly moisten the tube with a water soluble lubricant
- Insert the Provox LaryTube into the tracheostoma
- Attach the Provox LaryTube with Provox LaryClip, a Provox TubeHolder, or a Provox Adhesive
- Insert the Provox HME or FreeHands HME into the cassette holder



Provox LaryTube with Provox LaryClips

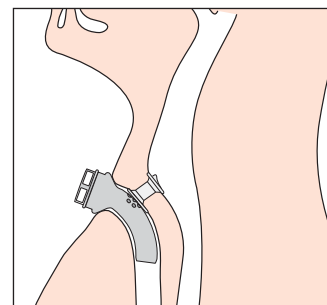

Hospital	_____
Department	_____
Clinician	_____
Contact information	_____



Provox® LaryTube™

When and How to Use

2. Cleaning the Provox LaryTube:


- Scrub thoroughly in a hot water-soap solution. Do not use synthetic detergents or oily based soaps. Rinse copiously in hot water and follow with a thorough rinse in distilled water.
- Clean the Provox LaryTube using a Provox TubeBrush. The Provox Brush can be used to clean the punctures of a Provox LaryTube with holes.
- Can be disinfected with boiling water for 5 – 10 minutes.
- Let the Provox LaryTube cool before reinsertion.

The Clinician will determine the appropriate size using the Provox LaryTube Sizer Kit.

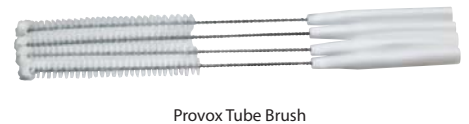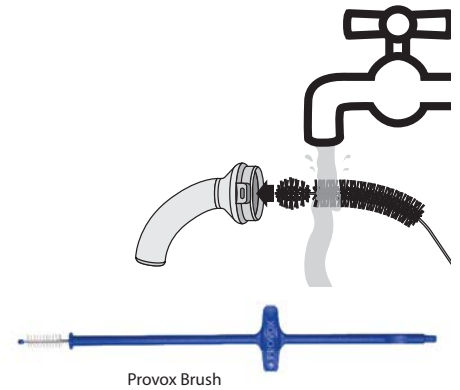

Provox®, Provox® HME, Provox® FreeHands HME® and Provox LaryClip™ are protected by patents, patents pending and design patents and design patents pending. Provox®, Provox® FreeHands HME® are registered trademarks owned by Atos Medical AB, Sweden. Provox® LaryTube™, Provox® LaryClip™ are trademarks owned by Atos Medical AB, Sweden.

This text does not replace nor does it set forth the complete contents of the User Manual and/or Prescriber Information for the products in this text, and is not a substitute for reviewing and understanding that important information. Therefore, before prescribing and/or using any of the products included in this text, please review the entire contents of the respective User Manual and/or Prescriber Information.
For more information contact us at: info.us@atosmedical.com • www.atosmedical.us • 1.800.217.0025

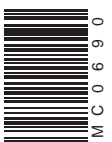

© Atos Medical AB, Sweden. 201107 MC0690-ThUS